

Electronic Balances

TX/TXB/TW Series







▶ 4200.00g

Clearly visible even in dark places

TX/TW_{series}

▶ 420.000g

Capable of performing delicate jobs





TX2202L, TX3202L, TX4202L

TX223L, TX323L, TX423L, TW223L, TW323L, TW423L

Wider applications ranging from smaller to larger objects



Value with no compromise





TXB_{series}

▶ 6200.0_g

▶ 620.00_g

Taking on weighing tasks everywhere

Reliable partner that creates weighing desires





TXB2201L, TXB4201L, TXB6201L, TXB6200L

TXB222L, TXB422L, TXB622L, TXB621L

Light, strong, high-quality design can be operated by batteries



Performance to allows quick weighing of smaller objects anywhere Creates weighing enjoyment





WindowsDirect Communication Function Only a cable is needed to send weighing data to a computer

TX / TW

Send balance data to Excel or other Windows applications without any data communication software installation required. All you have to add is one RS-232C cable. By combining standard AutoPrint functions with typical spreadsheet functions, even difficult applications can be easily automated. If you'd like to use "WindowsDirect" with "Windows 7" "Windows Vista", or USB port, please contact to our distributors





Can be used anywhere with

TXB

batteries.

Comparator function

TX/TW

Compare samples to target values or pass/fail criteria and clearly indicate the results.

RS-232C interface

TX / TW TXB

All models include standard RS-232C interface for easy integration with other devices and computers

Expanded Piece Counting function

Unit weights of up to 5 different samples can be easily entered, stored and recalled for use.



Streamlined design for anti-theft options

single touch

TXB

Industry-standard Kensington Security Slot accommodates the widest variety of optional anti-theft security devices.

Change among many weighing units and functions with a TX/TW

In addition to grams (g), weigh in %, number of pieces, ct, kg, mg, lb, oz, TTI, etc. or a custom conversion unit, more than 20 units in all. Change quickly from display of % or counting to gram weight (g) display.

Power saving function

TX/TW

When weighing operation ends, power automatically turns off after a fixed or pre-set time.

Fast and flexible calibration

TX/TW TXB

Balance display provides prompts for the calibration procedure and accepts your choice of external calibration weight.


Many useful functions included

TX / TW

- ·Alpha-numeric menu display
- •Up-front level indicator on display panel Simplified door cleaning
- ·Stainless steel pan
- ·Choice of weighing modes
- •Zero tracking function and settings ·Continuous data output function

The beginning of the new standard



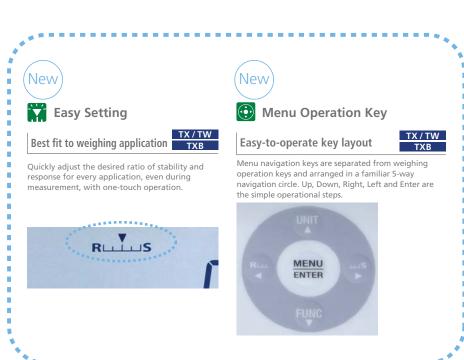



Touch-key Internal Calibration

TW

Press just two keys to calibrate the balance whenever calibration is necessary.
Calibration is very fast, taking only 15 sec.


Durable, high-performance aluminum alloy mass sensor

UniBloc Since 1989

TX/TW

The TX/TW series incorporates a one-piece aluminum alloy mass sensor technology (UniBloc), first introduced by Shimadzu for precision balances in 1989. It excels in performance, and resists deterioration and damage by ordinary impacts. The UniBloc's compact, uniform structure replaces 70 parts found in a conventional electromagnetic balance sensor assembly and ensures stable temperature characteristics, excellent response time and stable corner-load performance.

The UniBloc design permits a consistency of production that assures reliability and a long operational life.

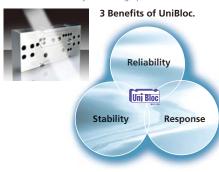


Fine measurements are assured with sliding glass door

windbreak system

TX/TW

TX/TW models with 1mg minimum display feature a sliding glass door type windbreak as standard configuration. The highly transparent system includes semi-transparent corners and maximum door opening capability.

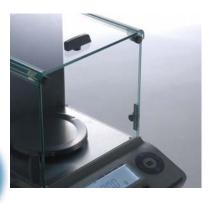

Illuminated display provides comfortable display visibility in all settings.

TX/TW TXB

A liquid crystal indicator with an attractive light blue backlight is standard on all models. The TXB series has an additional power-saving feature: the backlight automatically shuts off when idle and lights up again when weighing operations resume.



Specifications

TX/TW series Uni Bloc

17/ I VV SCITES CITES							
Model	TX223L / TW223L	TX323L / TW323L	TX423L / TW423L	TX2202L	TX3202L	TX4202L	
Capacity	220 g	320 g	420 g	2200 g	3200 g	4200 g	
Minimum Display	0.001 g	0.001 g	0.001 g	0.01 g	0.01 g	0.01 g	
Repeatability		0.001 g		0.01 g			
Linearity	±0.002 g						
Operating Ambient Temperature (°C)	5–40						
Temperature Coefficient Sensitivity (ppm/°C) (10–30°C)	±3						
Pan Size (mm) approx.	ø110 167(W) × 181(D)						
Main Body Dimensions (mm) approx.	206(W) × 291(D) × 241(H) 200(W) × 291(D) × 80(H)				(H)		
Weight (kg) approx.	Т	X: 3.8 kg TW: 4.2 k	g	2.8 kg			
Power Requirements			12V,	1A			

TXB series

Model	TXB222L	TXB422L	TXB622L	TXB2201L	TXB4201L	TXB6201L	TXB621L	TXB6200L
Capacity	220 g	420 g	620 g	2200 g	4200 g	6200 g	620 g	6200 g
Minimum Display	0.01 g			0.1 g				1 g
Repeatability	0.01 g			0.1 g				1 g
Linearity	±0.01 g	±0.02 g		±0.1 g	±0.2 g		±0.1 g	±1 g
Operating Ambient Temperature (°C)	5–40							
Temperature Coefficient Sensitivity (ppm/°C) (10–30°C)	±15 ±10 ±5 ±15 ±10 ±5		±	±20				
Pan Size (mm)	ø110				ø160	ø110	ø160	
Main Body Dimensions (mm) approx.	199(W) × 260(D) × 77(H)							
Weight (kg) approx.	1.5 kg							
Power Requirements	12V 1A or 6AA batteries							


TX/TXB/TW series options

•
Description
Printer EP-80
Printer EP-90
RS-232C cable
In-use protective cover
In-use protective cover for display



Extremely capable, but easy to operate

Dimensions

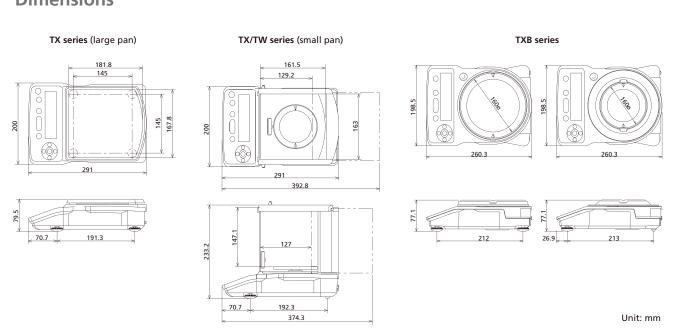



Shimadzu Corporation www.shimadzu.com/an/

Company names, product/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation or its affiliates, whether or not they are used with trademark symbol "TM" or "®".

Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services. Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

For Research Use Only. Not for use in diagnostic procedures.
The contents of this publication are provided to you "as is" without warranty of any kind, and are subject to change without notice. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.